

# Medcare 625 Ceiling Lift



Manufactured in the United States, Medcare Ceiling Lifts are the leading transfer equipment choice for patients with limited mobility or little weight bearing capacity. Medcare Ceiling Lifts are ideal for rooms where space is limited. Medcare's assortment of track options and components allow installation in a variety of room and building styles. Combined with Medcare's training program, ceiling-mounted lifts increase nurse and caregiver usage therefore increasing compliance with safe handling policies. With a weight capacity of 625 lbs., Medcare Ceiling Lifts dramatically reduce and eliminate transfer-related nurse injuries, and direct and indirect costs associated with these injuries.

#### Item Number:


323120 – Medcare Ceiling Lift 625 – lifts up to 625 lbs.

323102 – Medcare Ceiling Lift 450 – lifts up to 450 lbs.

#### Features:

- *Digital display time meter*
- *Manual or power traversing*
- *Return to charge option*
- *Continuous charge option or midtrack charger*
- *Visual/audible low battery indicator*
- *Emergency lowering system*
- *Built-in safety stops*
- *Pneumatic hand control*
- *Scale option*
- *Automatic battery conservation*
- *Soft start/stop for comfort*
- *Heavy duty batteries for more lifts*
- *Efficient battery charge cycle*
- *2 year warranty (excluding batteries)*

“Medcare Ceiling Lifts provide a safe, smooth and effortless transfer for patients and caregivers.”

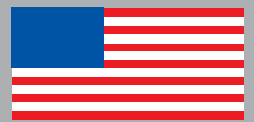


MADE IN THE USA



# Medcare 625 & 450 Ceiling Lifts

Model #323120 – lifts up to 625 lbs. Model #323102 – lifts up to 450 lbs.

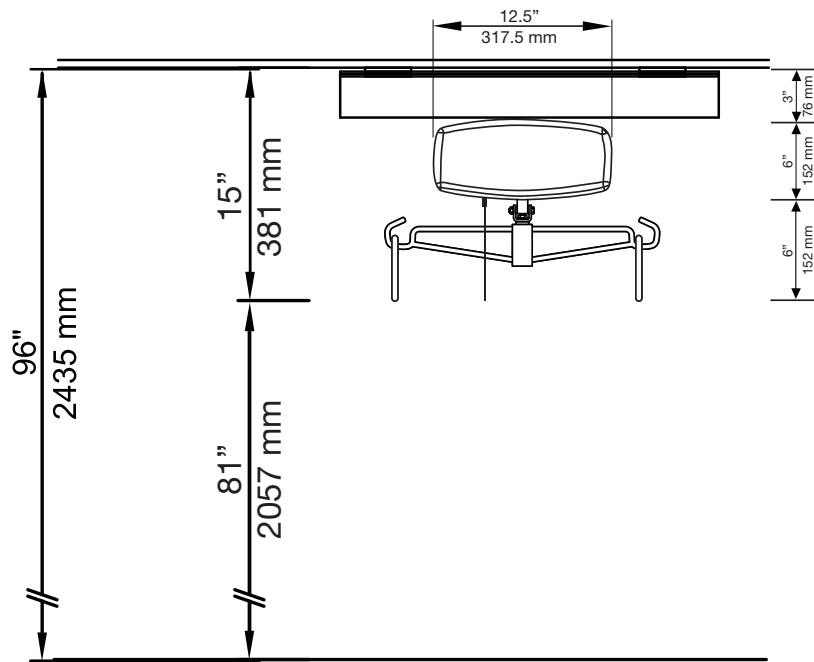

## Technical Specifications

<b>Lift Motor:</b>	24 VDC (All models)	<b>Lift Case:</b>	Flame Retardant ABS
<b>Traverse Motor:</b>	24 VDC (Optional at time of purchase)	<b>Hand Control:</b>	Pneumatic
<b>XY Traverse Motors:</b>	24 VDC (Optional at time of purchase)	<b>Lifting Height:</b>	Up to 96" (Motor strap is 81" long when fully extended)
<b>Charger Input:</b>	120 VAC, 1.0 Amps	<b>Hoist Weight:</b>	21-23.5 lbs.
<b>Charger Output:</b>	24 VDC, 1.5 Amps	<b>Maximum Load:</b>	625 lbs.
<b>Batteries:</b>	24 VDC (2 x 12 VDC), 5.0 AH, Sealed Lead Acid	<b>Warranty:</b>	2 Years (excludes batteries)

Med  
HANDLE WITH CARE.

### Hanger Bar Stopper Option

The Medcare Hanger Bar Stopper prevents sling loops from sliding off of the hanger bar hooks.



**Item Number:**  
400611

## PRODUCT HIGHLIGHTS

- Digital display running time meter along with audio and visual battery life indicator
- 96" lifting range allows lifting from the floor and over high thresholds
- The slack strap sensor ensures the lift strap will not continue to unwind if weight distribution is not even.
- Manual lowering in the absence of power
- 4-point, non-rigid hanger bar design allows sling to conform to patient for added comfort



### Medcare Products

151 East Cliff Road  
Burnsville, MN 55337

Toll-free: 800-695-4479

Fax: 952-894-7153

www.medcarelifts.com

© 2009 Medcare Products Inc.  
Specifications are subject to change without notice. All products referred herein are registered trademarks or trademarks of their respective holders.